


New Agro Resea Noorddijk 3 2391 CE Hazers The Netherlands T. +31 (0)172 419 F. +31 (0)172 419	woude dorp 9 449	PEAK RESULTS	info@newagroresearch.nl www.newagroresearch.nl IBAN NL 64 ABNA 0563 8714 66 VAT nr. NL8563.46.780.B.01 Chamber of Commerce 65990692
Custome	r Details	Sample details	
Company Contact Address Zip & City Telephone Email	Procare Nutraceuticals BV Martyna Kriukovska Noorddijk 3 2391 CE, HAZERSWOUDE DORP 0172 419449 info@procare.nl	NAR Barcode Sample name Type of sample Customer Batch Number Customer intermal code Sample taken by representative for whole batch comment	20180830 - 04 PN - Hempseed oil 10% Oil 20180830 - 04 / 0051 New Agro Research Yes
Injection Deta			
Instrument	Agilent1100	Sample Weight in Mg	78.3000
Method	201808 IM 420	Dilution Factor:	20.0000
Sequence Name	201808230	Injection Volume:	10.00
Injection Name:	PN - Hempseed oil 10%	Injection Date/Time:	31-aug-18 16:34
Injection Number: Vial Number:	6 Vial:4	Run Time (min):	53.99

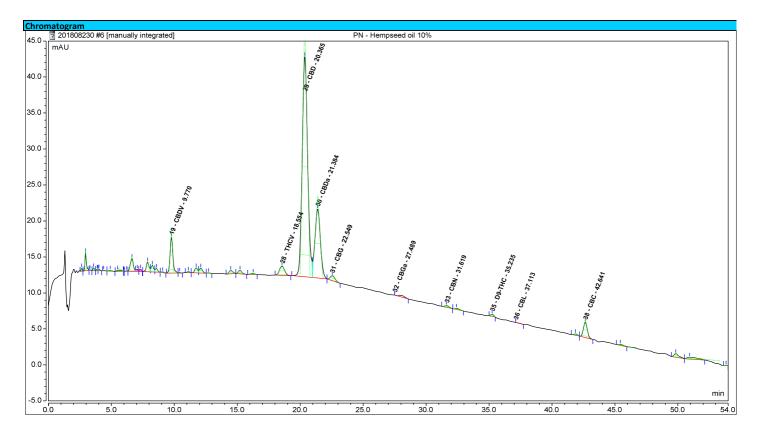

No.	Full Name	Peak Name	Retention Time	Amount	Amount	Amount
			min	%	mg/G	ррт
19	n.a.	CBDV	9.77	0.78	7.8	7778
28	n.a.	THCV	18.55	0.36	3.6	3604
29	n.a.	CBD	20.36	9.25	92.5	92471
30	n.a.	CBDa	21.38	1.59	15.9	15919
31	n.a.	CBG	22.55	0.19	1.9	1921
32	n.a.	CBGa	27.49	0.04	0.4	392
33	n.a.	CBN	31.62	0.03	0.3	301
35	n.a.	D9-THC	35.24	0.05	0.5	484
36	n.a.	CBL	37.11	0.00	0.0	12
38	n.a.	CBC	42.64	0.20	2.0	2035

Total THC's		0.05	%		
Ratio	THC : CBD		220		


Following cannabinoids are not present in sample:

D8-THC Should have appeared at 36.263 THCa Should have appeared at 44.063

Acidic cannabinoids do not decarboxylate in to neutrals completely and will not have the same effect as their neutral forms, CBD and THC. That is why we use a correction factor of 0.877 for the total percentages above in the blue bar. If you add the percentages of CBD and CBDA in table 1, this will be higher than the total percentages.

Total of detected Cannabinoids

Relative percentage of the components over the total of the detected cannabinoids. (these are higher than the percentages in the analyzed sample, because all other substances like plant material, proteins, fats etc) The real value is in the table "Identified Peak Results" above

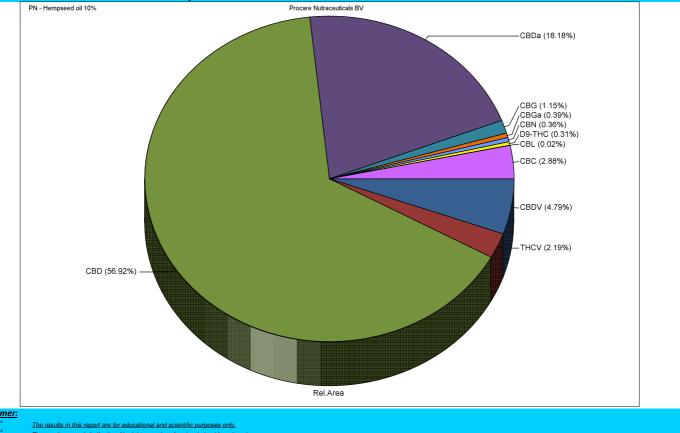

<u>Ine executive laboratory and its personnel accept no liability whatsoever for errors in the reportin</u> <u>nor for material and/or immaterial damage resulting from the services it provides.</u> Instrument: Agilent1100 Sequence: 201808230

20180725-1 Extern(1)/Peak Analysis

All rights reserved New Agro Research BV